

# Astor Compact DCII

Dementia Care Height Adjustable Bath with Powered Seat

### Highlights:

- **Dementia Care** features and lockout
- **Height adjustable** to minimise carer handling strain
- **Integrated powered seat system** for safe transfer (detachable option)
- **Seat transfer trolley** for easy bathroom to bedroom transfers (with detachable seat option)
- **Robust** panels and leg covers
- **63" and 67"** options




The **Astor Compact DCII Bath** is the first bath specifically designed for use within a dementia care setting. The unique auto lockout ensures that the water cannot be accidentally turned on and that the bath features are immobilised. The bath has a built in auto-fill system to enable the carer to set the bathwater running unattended, without the risk of overflow.

The unique shape of the Compact bath range and height adjustable design give maximum comfort to the bather and minimal strain and ease of use for the carer. With a proprietary bearing system, robust panels, and high grade 304 stainless steel parts, the Compact range has been built to the highest standard, **designed for longevity**.

The **integrated powered seat**, with its clip and lock lap strap, nursing arms and obstacle sensor, enables bathers to transfer safely and comfortably. On the Compact DCII, there is an option for the seat to be of a contrasting colour to help bathers with perception issues to sit down more assuredly. The seat is shaped to give optimal carer access, especially with the commode seat option. There is also a detachable seat option where the seat can be completely removed, and a **headrest fitted** to enable comfortable bathing when the seat is not in use. A transfer trolley is available for safe and smooth transfers.

The bath water controls are fitted with a TMV3 temperature control with a digital display showing the temperature of the water leaving the spout.

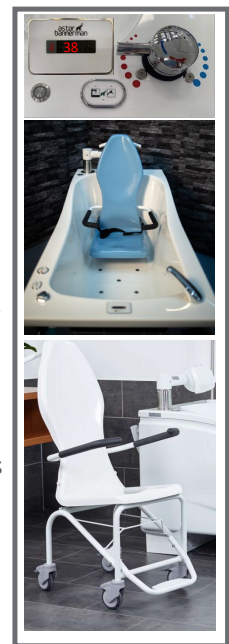

### Standard Features:

- Maximum Working Load (MWL) 24 st / 336 lb
- Powered seat with lap strap, nursing arms and obstacle sensor
- Bath filler system with built in TMV3 temperature regulation
- Autofill sensor to prevent overflowing
- Emergency battery backup
- Length Options: 63" or 67" (x 30")
- Working height range: 31" - 39"
- 3 Year warranty (*Terms and conditions apply*)

### Upgrade Features:

- Coloured seat
- Air spa system
- Bluetooth sound system
- Colour changing LED light system
- Detachable seat system with head rest
- Seat transfer trolley
- Commode seat

# Product Specification Sheet: Astor Compact DCII




<b>Model Name:</b>	<b>Astor Compact DCII</b>
<b>Product Description</b>	Dementia Care height adjustable bath with powered seat
<b>Height Range</b>	31" - 39"
<b>Key Dimensions</b>	63" or 67" length x 30" width

## Features:



## Options:

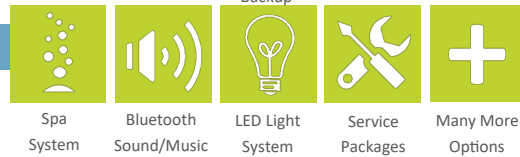

## Application:



## Technical Drawings:

All dimensions in inches:

