

Astor CT Nivano Wall Mounted Changing and Showering Table



Highlights:

- Height adjustable to give safe and comfortable carer working height
- Suitable for use in shower areas and wet rooms
- Durable non-rip high tension vinyl stretcher - comfortable and easy to clean
- 150 kg SWL (220 kg option)

A popular choice among families and Occupational Therapists for many years, the CT Nivano is comfortable, reliable and looks sleek. It is simple to operate from a handset control and comes in a range of sizes from 1000 mm up to 1900 mm to suit both paediatric and adult applications.

The durable stretcher material is comfortable as well as quick and easy to clean, saving precious time. This, along with part stainless steel construction gives many years of worry free use, even in wet and humid or multi-user environments, especially when combined with ongoing maintenance for extra peace of mind.

The 600 mm travel range enables the user to be gently lifted to a safe and comfortable working height for washing, showering, changing and dressing. The platform stows neatly when not in use, held in place with twin gas struts, and features battery back-up for safe lowering in the event of a power failure.

The CT Nivano can be tailored for you with the addition of optional extras adding comfort, convenience and unique customisation.



Standard Features:

- SWL 150 kg/24 st
- Height range 300mm - 900mm (600mm)
- Height adjustable (mains electric)
- Handset operation
- Battery backup
- Floor and ceiling hoist access
- Part stainless steel construction
- Durable non-rip high tension vinyl stretcher
- Protrudes just 242 mm from wall when stowed
- 3 Year warranty (*Terms and conditions apply*)

Optional Extras:

- Stretcher sizes from 1000-1900 x 700 mm
- SWL 220 kg/35 st (1700-1900 models)
- Battery operated model (CTB2)
- Fold away safety guard
- Protective padded guard bumper
- Anti-splash skirt
- Tilting head/backrest (1700-1900 models)
- Protective foam wedge support
- Removable mattress
- Special fully coloured models available
- Range of service & maintenance packages:

Product Specification Sheet:

Astor CT Nivano Changing & Showering Table



Model Name:	Astor CT Nivano
Product Description	Wall mounted changing and showering table
Product Codes	150 kg SWL: 00.038.11.1000 - 00.038.11.1900 (1000-1900 mm) 220 kg SWL: 00.038.12.1700 - 00.038.12.1900 (1700-1900 mm)
Mattress Size	150 kg SWL: 1000 mm - 1900 mm x 700 mm 220 kg SWL: 1700 mm - 1900 mm x 700 mm
Height Range	300 mm - 900 mm x 700 mm (600 mm)

Features:

Height Adjustable	Mains Electric	Battery Backup	150kg/24st SWL	3 Year Warranty

Options:


Battery Model	1000 - 1900mm	220kg/35st SWL Model	Tilting Head/Backrest	Service Packages	Many More Options

Application:

Private Home	School	Leisure Centre	Care Home	Nursing Home	Respite Home	Hospital	Changing Places
--------------	--------	----------------	-----------	--------------	--------------	----------	-----------------

Technical Drawings:

All dimensions in mm:

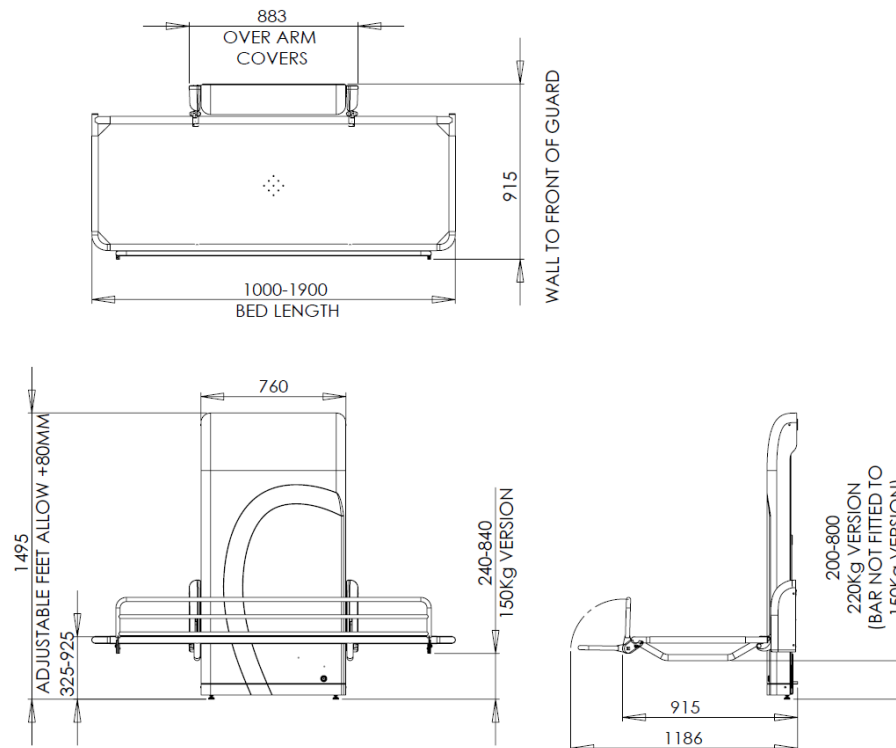

Astor-Bannerman, Andoversford Business Park, Cheltenham GL54 4LB
 sales@astorbannerman.co.uk **01242 820820** www.astorbannerman.co.uk