


Product Specification Sheet:



Astor CTX2-CP

Wall Mounted Changing and Showering Table for Multi-User Environments

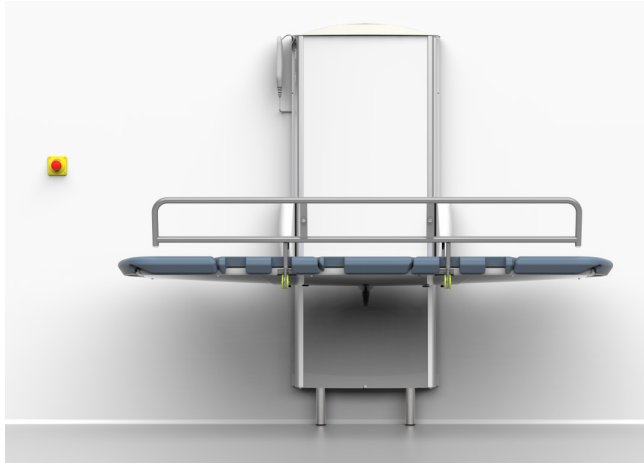

Highlights:

- 200 kg / 31½ st maximum working load (MWL)
- Discreet integrated water collection tray
- Calming ambience through integral lighting
- Comfortable polyurethane padded mattress
- Flexible support with tilting headrest
- Guards, bumpers and wall mattress available


The Astor CTX2-CP is the ultimate wall mounted changing table combining many practical and aesthetic features to allow the complex needs of daily care giving to be conducted in a comfortable and relaxed manner.

The height of the bed is easily adjusted from a low of 300 mm/11¹³/₁₆" to 1040 mm/40¹⁵/₁₆" via the handset to minimise carer strain, and to optimise manoeuvrability. The handset also controls the integrated lighting which gives a calm and relaxed environment for transferring and changing. Tilting head rests give enhanced comfort, combining with other features to enhance the personal care routine.

The Astor CTX2-CP includes a number of added safety features including an emergency stop button, battery back up in case of a power failure, and a power status light.

Our changing table has been designed to withstand high usage in wet and humid environments with the inclusion of stainless steel parts, a durable and comfortable moulded polyurethane bed with drain holes and an integrated water collection tray, making it a suitable addition to wetrooms and ensuites.

The mattress stows neatly when not in use to maximise room space, held in place with twin gas struts.

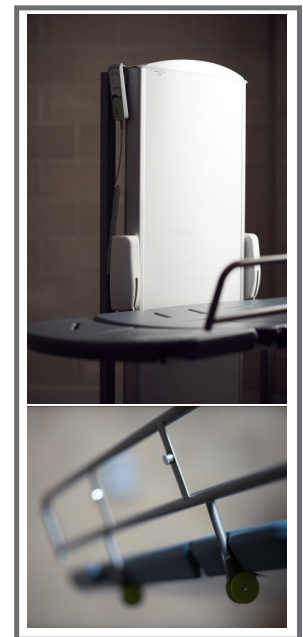

Standard Features:

- MWL 200 kg / 31½ st
- Height adjustable range:
300 mm - 1040 mm (740 mm) / (11¹³/₁₆" - 40¹⁵/₁₆")
- Handset operation (Mains electric: Battery backup)
- Rust free stainless steel bed frame & components
- Moulded polyurethane padded bed/platform
- Integrated water collection tray with drain hose
- Emergency stop button
- Battery back up and power status light
- Integrated LED light
- Integrated tilting headrest
- 3 Year warranty (*Terms and conditions apply*)

Optional Extras:

- Bed lengths 1500 or 1800 mm
59¹/₁₆" or 70⁷/₈"
- Battery operated model
- Fold away front safety guard
- Twin position end safety guards
- Protective padded guard bumper
- Removable mattress
- Wall cushion
- Range of service & maintenance packages




Astor CTX2-CP Changing and Showering

Model Name:	Astor CTX2-CP
Product Description	Wall mounted changing and showering table with front guard
Product Codes	00.036.17.1528 (1500 mm / 59 ¹ / ₁₆ ") 00.036.18.1828 (1800 mm / 70 ⁷ / ₈ ")
Mattress Size	1500 mm - 1800 mm x 700 mm / 59 ¹ / ₁₆ " or 70 ⁷ / ₈ " x 27 ¹³ / ₁₆ "
Height Range	300 mm - 1040 mm (740 mm) / 11 ¹³ / ₁₆ " - 40 ¹⁵ / ₁₆ "

Features:

							
Height Adjustable	Mains Electric	Battery Backup	200kg/31½st MWL	Tilting Headrest	LED Strip Light	3 Year Warranty	Front Guard

Options:

				
Battery Model	1500mm or 1800mm	Removable Mattress	End Guards	Service Packages

Application:

Private Home	School	Leisure Centre	Changing Places 	Respite Home	Hospital	Care Home
--------------	--------	----------------	---	--------------	----------	-----------

Technical Drawings:

All dimensions in mm:

